



NECK MASSAGER

User Manual Model: MG11



www.iyesido.com

Guangzhou Guan Su Ge Technology Co., Ltd.
Room A06, Floor 2, T28 Class, No. 16, East Of Hao Ke Zhou
Street, Shixi Village, Hai Zhu District, Guangzhou, China.

Thank You For Purchasing MG11 Product!
The Brief Description Of The Product Assures You
To Operate Freely, Fully Embodies The Value Of
This Product, Please Read Carefully Before Use.




List Of Package

1. Neck Massager
2. User Manual
3. Charging Cable





Function Diagram

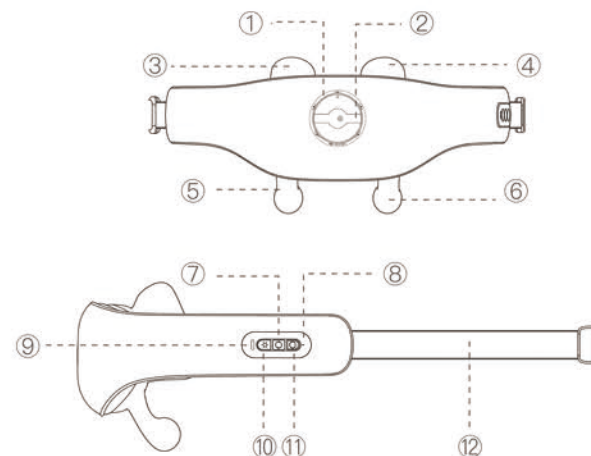

- | | |
|--------------------|--------------------------|
| 1. Neck Distance | 7. Mode |
| 2. Adjustment Knob | 8. Indicator Light |
| 3. Massage Head | 9. Charging Port |
| 4. Massage Head | 10. The Heating Function |
| 5. Massage Head | 11. Power Switch |
| 6. Massage Head | 12. Shoulder Strap |

How To Use 1

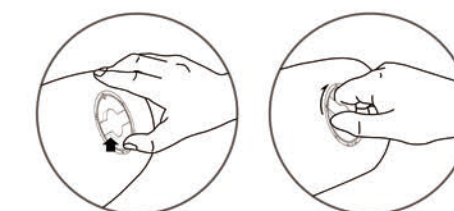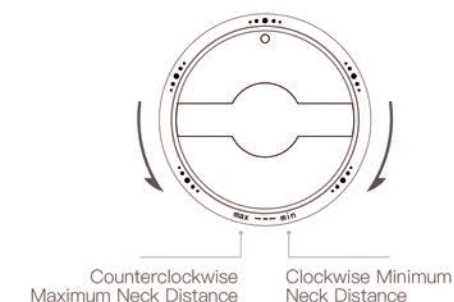

Grip The Knob And Turn It In The Indicated
Direction To Adjust The Neck Distance.

How To Use 2

Correct Usage Instructions

Preparation Before Use
Please Check If The Product And Accessories Are Complete And Undamaged.

- **Power On**
Press And Hold The Power Button For 2 Seconds To Start The Product.
- **Mode Selection**
The Default Mode On Startup Is Mode 1. Press The Mode Button Briefly To Switch Between Mode 1, Mode 2, And Mode 3 In A Cycle.
- **Heating Mode**
The Default Heating Level On Startup Is High. Press The Heating Button Briefly To Switch Between High And Low Settings. Press And Hold The Heating Button For 2 Seconds To Turn Off The Heating.
- **Neck Distance Adjustment**
The Product Defaults To Medium Neck Distance When Shipped. Turn Clockwise To Decrease The Neck Distance And Counterclockwise To Increase It.
- **Power Off**
Automatic Power Off: The Product Will Automatically Turn Off After 10 Minutes Of Use. Manual Power Off: Press And Hold The Power Button For 2 Seconds To Turn Off The Product.


How To Use 3

Charging Instructions

1. The Product Cannot Be Used While Charging.
2. Insert The Charging Cable Plug Into The Charging Port On The Product.
3. When The Product Is In Operation And You Connect It To The Power Source For Charging, The Product Will Activate Its Self-protection Mode And Automatically Stop Working.
4. When The Battery Level Drops Below 20%, The Indicator Light Will Flash Red.
5. After Charging Is Complete, The Indicator Light Will Remain Solid White.
6. It Is Normal For The Main Unit And Power Adapter To Become Slightly Warm During Charging.


End Of Charging

After Charging Is Complete, Please Unplug The Power Source.

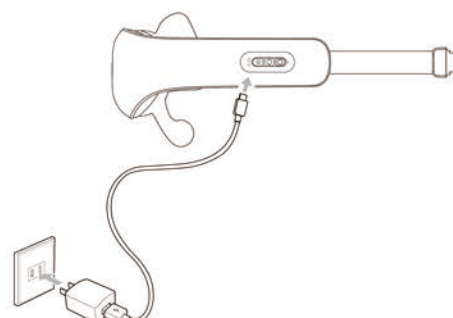

Product Parameters

Wearing Method

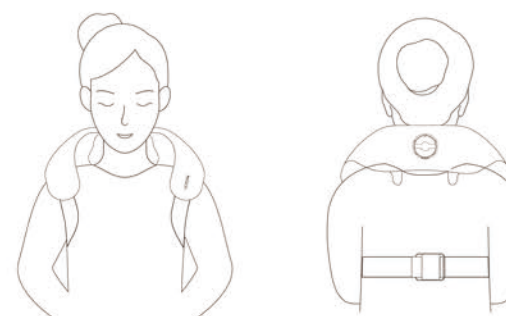

Model: MG11
Material: PU+ABS
Rated Power: $\leq 10W$
Battery Capacity: 2400mAh
Charging Time: About 3 Hours
Usage Time: 10 Min Per Cycle

Charging Interface: Type-C
Heating Temperature: $41^{\circ}C-45^{\circ}C$
Continuous Use: About 150 Min
Rated Voltage: DC 7.4V
Charging Current: DC 5V/1.5A
Noise: 50dB

Safety Instructions

Storage

If The Product Is Stored For More Than 3 Months, Open The Box, Inspect The Machine, And Recharge It Before Use. When Storing, Avoid Placing The Product Directly On The Floor, Against Walls, Or Near The Ceiling. Ensure Good Ventilation Indoors And Avoid Exposure To Strong Sunlight And Corrosive Gases.

Transportation

The Product Packaging Includes Simple Shockproof Features And Is Suitable For Air, Rail, Road, And Sea Transport. Avoid Exposure To Rain, Snow, And Splashes, As Well As Inversion And Collision.

Power Check

Check The Power Supply; Do Not Operate The Product If The Power Voltage Is Outside The Specified Range.

Maintenance

Take Measures To Prevent Moisture, Rodents, And Pollution. Clean The Product While It Is Turned Off. After Use, Place The Device In A Safe Location And Avoid Placing Heavy Objects On It. If A Fault Occurs, Seek Repair From A Qualified Technician; Do Not Disassemble The Product Yourself.

Disposal And Recycling

This Product Should Not Be Disposed Of With Household Waste. When You Decide To Dispose Of The Product, Please Follow Local Environmental Laws And Guidelines For Disposal.