



BACK MASSAGER

User Manual

Model: MG15




www.iyesido.com

Guangzhou Guan Su Ge Technology Co., Ltd.
Room A06,Floor 2,T28 Class,No.16,East Of Hao Ke Zhou
Street, Shixi Village, HaiZhu District,Guangzhou,China.

Product Display

The Brief Description Of The Product Assures You
To Operate Freely, Fully Embodies The Value Of
This Product, Please Read Carefully Before Use.

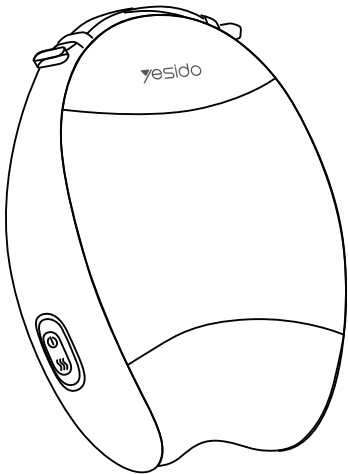


Thank You For Purchasing MG15 Product!

Product Parameters

Model: MG15
Input: 5V/2A
Material: PU+ABS+Fabric+Cotton
Noise: <60dB
Rated Voltage: 7.4V
Rated Power: 32W
Battery Capacity: 2000mAh
Charging Time: About 3 Hours
Continuous Use: About 90 Min
Usage Time: 10 Min Per Cycle
Charging Interface: Type-C
Heating Temperature: 30°C~40°C
Massage Head Speed: 26~30 RPM


Function Introduction

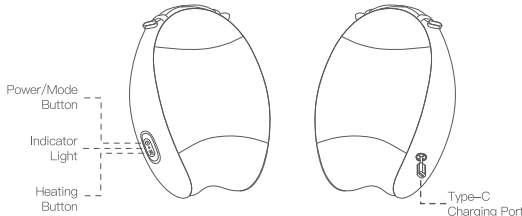


- 1.Power Button**
Long Press “” To Power On, The Green Light Will Keep On And Enter Mode 1;
Long Press “” To Power Off.
- 2. Mode Button**
Start Up To Default Mode 1, Press “” To Switch To Mode 2.
M1: Kneading + Tapping(10 Min)
M2: Tapping(10 Min)
- 3. Heat Button**
The Default Mode Does Not Heat When Start Up. Press “” To Turn The Heating On Or Off. When The Heating Key Is Switched On, The Green Indicator Blinks Twice. Massage Head Position Bright Red Light. The Heat Is Off And The Red Light Is Off At The Same Time.
- 4. Indicator Light**
Working, The Green Light Is On For A Long Time And Blinks Green Twice.
Low Power, The Red Light Blinks Until The Power Is Turned Off.
Fully Charged: The Green Light Is On. (Can Not Be Turned On When Charging)
- 5. Type-c Charging Port**
Connect The Type-C Cable To The Charging Port To Charge.

Product Introduction

Mode:

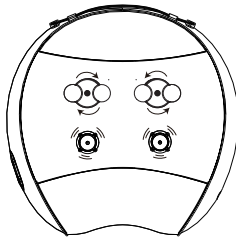


Mode 1
Knead + Tapping

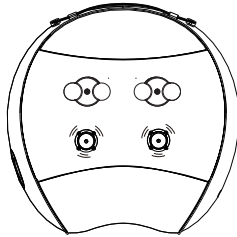

Mode 2
Tapping mode

Exploded Drawing Of Battery

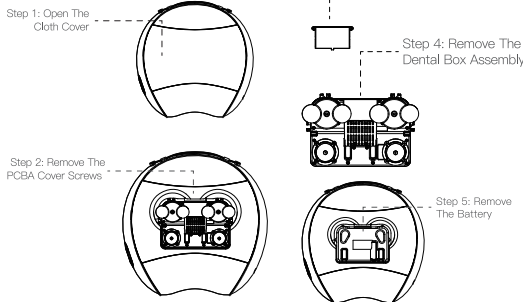

Precautions

Persons Not Applicable

Never Use It With The Following Medical Electronic Instruments

1. Implantable Medical Electronic Instruments Such As Cardiac Pacemaker.
2. Life-sustaining Medical Electronic Devices Such As Artificial Heart-lung Machine.
3. Forbidden For Pregnant Women.

Precautions

If You Feel Unwell During Or After Use, Discontinue Use And Consult A Doctor;This Product Contains Lithium Batteries, Please Do Not Disassemble, Modify, Immerse The Product In Water Or Throw It Into Fire, Otherwise It May Cause Damage To The Product And Explosion; Please Consult A Physician Before Using This Product. Otherwise, Accidents Or Physical Discomfort May Occur

1. Patients With Abnormal Heart And Lumbar Nerves.
2. Patients With Severe Hyperosteoegeny.
3. People With Abnormal Spine Or Back Curvature Caused By Accidents Or Diseases.
4. People With Severe Blood Circulation Disorder At The Waist.
5. Patients Who Cannot State Their Wishes.
6. Pay Attention To Use This Device If You Are Insensitive In Heating. Please Do Not Use It For Any Purpose Other Than Massage Recorded In The Instruction Manual,Otherwise, Accidents And Faults May Occur. Please Do Not Disassemble, Repair Or Modify This Product Without Authorization, Otherwise, Faults Or Electric Shock Accidents May Occur.

Use Environment

1. Please Do Not Use This Machine In Damp Or Dusty Environment, So As To Avoid Electric Shock Or Mechanical Faults.
2. Please Do Not Place This Machine In Direct Sunlight, And Near Electric Heaters Or Stoves.
3. Please Do Not Use This Machine At Room Temperature Above 40°C

Precaution & Maintenance

If You Encounter Problems When Using The Product, Please Refer To The Following Tips And SolutionsIf The Product Is Still In A Malfunctioning State Or Cannot Be Used, Please Turn Off The Power. Usersshould Not Disassemble Or Repair The Product By Themselves. Damage To The Product Caused Bydisassembly Or Repair Will Invalidate The Product Warranty, Our Company Will Not Be Responsible Forany Damage Arising From This Situation.

Conditions	Possible Reasons/treatment Methods
Unable To Start Massager	Please Confirm Whether The Battery Is Low; Please Confirm That The On-off Button Has Been Pressed.
Unable To Charge Or No Charge Lamp Display	Please Check Whether The Type-c Charging Cable Is Connected To The Charging Port Of This Machine.
Massager Suddenly Stop During Work	Check If The Battery Is Low. Check If 10-minutes Timer Have Reached. Please Don't Continue To Use It, And Contact The Maintenance Service Center If The Problem Still Exists.

Precaution

Danger

The Battery Used In This Product Is A Special Rechargeable Battery, Which Cannot Be Used In Other Products. Do Not Charge The Rechargeable Battery That Has Been Taken Out;
Do Not Throw Into Fre Or Make It Overheated.
Do Not Pierce With Nails, Impact With Gravity Or Dismantle,
Do Not Connect The Positive And Negative Electrodes Of The Battery Directly With Metal;
Please Do Not Store It With Necklaces, Hairpins, Etc. Please Don't Charge. Use Or Place It Under High Temperature Places Near Freside Or Under Hot Sunlight. Or Place It In High-temperature Places Such As Freside Or Hot Sun. The Above Conditions May Cause Heating, Fre, Rupture And Other Phenomena;

Warning

Please Keep The Removed Rechargeable Battery Out Of The Reach Of Infants To Prevent Accidental Swallowing; If Is Swallowed By Mistake, Please Seek Medical Attention Immediately;when The Battery Fluid Leaks, Do Not Touch The Fluid With Hands Without Protection.and Handle It As Follows:
When The Battery Fluid Gets Into The Eyes, There Is A Danger Of Blindness. Do Not Rub Your Eyes Directly With Your Hands. Wash Them With Clean Water Immediately And Seek Medical Attention;
When Battery Fluid Gets On Your Body Or Clothes, It May Cause Skin Irritation Or Injury. Please Rinse With Clean Water. If You Still Feel Unwell, Seek Medical Attention Immediately;
Only The Product Can Be Decomposed When It Is Discarded. Failure To Comply May Result In Fire, Electric Shock, Injury And Other Accidents;

Warranty Terms

1. This Product Comes With A One-year Free Warranty From The Date Of Purchase. Please Present This Warranty Card When Requesting Repairs.
 2. If The Warranty Period Has Expired, Repairs Can Be Made At The Company's Designated After-sales Service Center For A Fee.
 3. Damage Caused By Human Error Or Appearance Damage Is Not Covered Under The Warranty.
- High-quality After-sales Service Is Provided To Protect Your Legal Rights And Interests. To Facilitate Repair Services, Please Keep The Warranty Card In A Safe Place. The Company May Charge For Parts Based On The Actual Situation During Repairs.

Product Warranty Card

Product Warranty Card

Purchase Date	Year	Month	Day
Product Name			
User Address			
Name			
Phone Number			
Dealer			
Failure Cause			

The Final Interpretation Of The Product Belongs
To Guan Su Ge Technology.