



HEATED EYE MASK

User Manual Model: MG13



www.iyesido.com


Guangzhou Guan Su Ge Technology Co., Ltd.
Room A06, Floor 2, T28 Class, No. 16, East Of Hao Ke Zhou
Street, Shixi Village, HaiZhu District, Guangzhou, China.

Thank You For Purchasing MG13 Product!
The Brief Description Of The Product Assures You
To Operate Freely, Fully Embodies The Value Of
This Product, Please Read Carefully Before Use.



*Display Of Appearance

Function Diagram

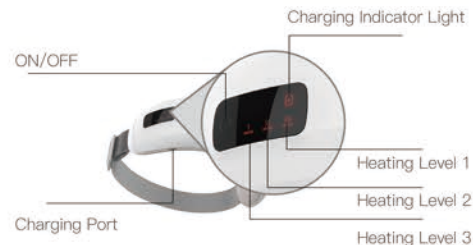

List Of Package

Heated Eye Mask
User Manual
Charging Cable



Product Parameters

Power On/Off: Press And Hold For 2 Seconds.
Switch Mode: Press The Power Button Once After Turning On.



Mode Function
Heat Therapy Mode:
Heat Level 1 42°C $\pm$ 2°C
Heat Level 2 44°C $\pm$ 2°C
Heat Level 3 46°C $\pm$ 2°C

Model: MG13
Material: ABS
Rated Voltage: DC 5V
Rated Power: $\leq$ 2W
Charging Current: 5V $\approx$ 500mA
Charging Time: Type-C About 2 Hours
Battery Capacity: 3.7V 600 mAh
Working Timer: 15 Minutes Per Use
Continuous Use: About 120 Minutes
Mode: 42°C $\pm$ 2°C, 44°C $\pm$ 2°C, 46°C $\pm$ 2°C (Perceived),
3-level Heat Therapy Mode

Note:

1. Connect The External Power Cable To Charge. The Indicator Light On The Device Will Turn On, And It Will Turn Off When Fully Charged. The Machine Cannot Be Used While Charging.
2. All Modes Have Fixed Functions, And The Timer And Features Are Not Adjustable.
3. If Not Used For A Long Time Or If The Battery Is Low, Please Charge At Least Once A Month To Avoid Damaging The Battery And Extending Its Lifespan.

Safety Instructions

- Do Not Use The Eye Massager While Driving.
- Do Not Use The Eye Massager In Excessively High Temperatures Or Humidity.
- Do Not Use If You Have Symptoms Of Glaucoma, Cataracts, Or Retinal Detachment.
- Use With Caution If You Have Had Eye Surgery; Consult A Professional Physician. If You Feel Uncomfortable During Use, Please Stop Immediately.
- Use With Caution If You Have Had A Stroke, Hypertension, Or Hypotension; Consult A Professional Physician. If You Feel Uncomfortable During Use, Please Stop Immediately.
- Those With Mobility Issues, Physical Or Mental Health Conditions, Or Cognitive Impairments Should Use With Supervision. If You Feel Uncomfortable During Use, Please Stop Immediately.
- Avoid Using The Eye Massager After Consuming Alcohol.
- Ensure No Contact Lenses Are In Your Eyes Before Using The Eye Massager.
- It Is Recommended To Wear The Massager On Your Eyes Before Turning It On To Prevent Overinflation.
- This Product Is Not Intended For Diagnosis, Treatment, Or Any Commercial Use.
- For Cleaning, Use A Damp Cloth With Neutral Cleaning Solution Or Disinfectant; Avoid Cleaning With Plain Water.
- This Product Has A Heating Function; Please Be Cautious If You Are Sensitive To Heat.